



$Gaius\ DUO, installation\ instructions.$


17782 wall drain 17783 floor drain

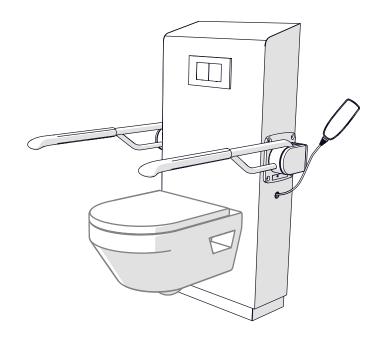


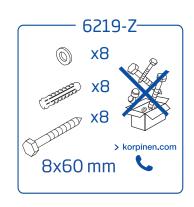

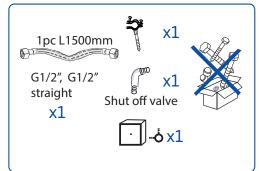

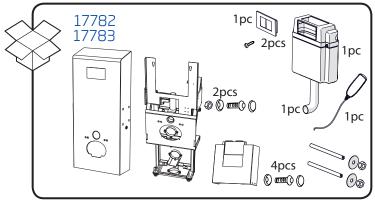

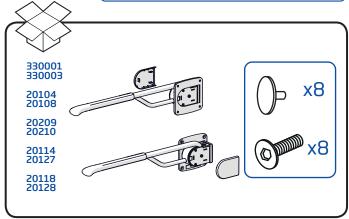

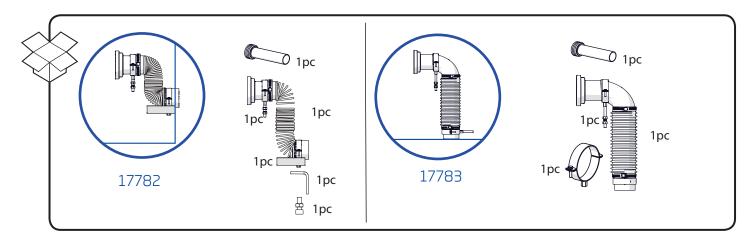



Gaius DUO,


Attachment points, 8 pcs Seinäkiinnitys 8 kpl

17782, 17783 17785, 17786

Gaius DUO,

Water and electric Vesi ja sähköpisteet

17782, 17783

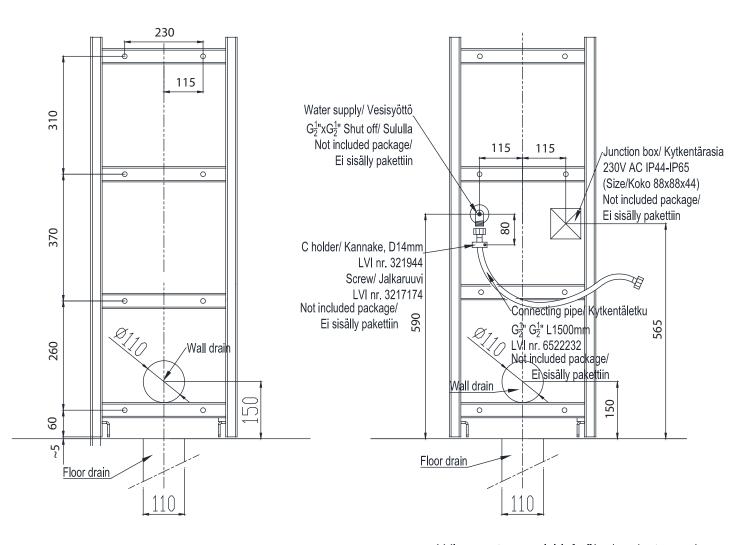

LVI nr. go to www.lvi-info.fi/en/product-search and find LVI number "Search content" LVI nr. mene www.lvi-info/tuotehaku ja etsi LVI-numerolla tuote

Structure

The product has undergone static load testing at 500 kg. The size of the wall mounting fixtures is 700x1150 mm. The fixtures must be affixed to solid structures. The fixtures must be mounted in accordance with the measurements specified in the diagram.

Electrics

The electrical junction box must be positioned to the right of the mid-line, as shown in the drawing. IP44 or higer junction box (not included)

The base comes with a 230V AC 10A socket.


A residual-current device is recommended for the 230 V AC power supply.

A qualified electrician must perform the connection.

HVAC

The water connection must be positioned to the left of the mid-line, as shown in the drawing.

- 1. Mid-line
- 2. Wall
- 3. Electrical connection
- 4. Water connection. Water supply with a shut-off valve equipped with a G1/2"external thread. Thread downwards
- 5. Wall mounting fixtures
- 6. Finished floor
- 7. Wastewater connection on the wall or on the floor (Ø110 mm). The pipe is cut at the level of the wall.
- 8. Reinforce area, plywood 36mm or over, solid wood 45mm or over.

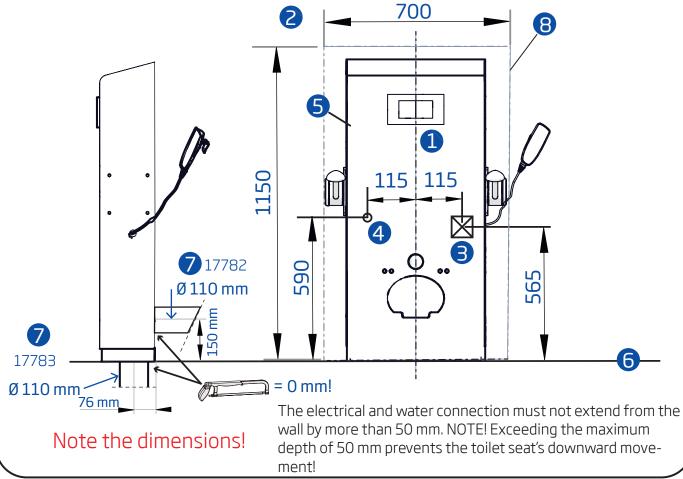



Notice!


General instructions:


Use motor up and down to help installation.


If frame don't go down by motor, press moving frame by hand.

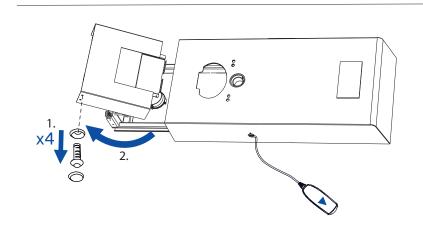


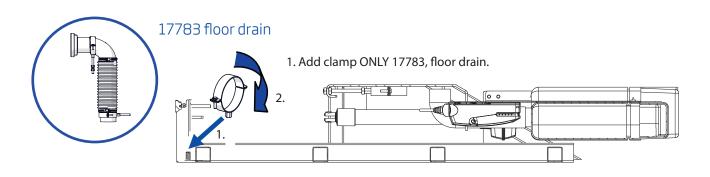


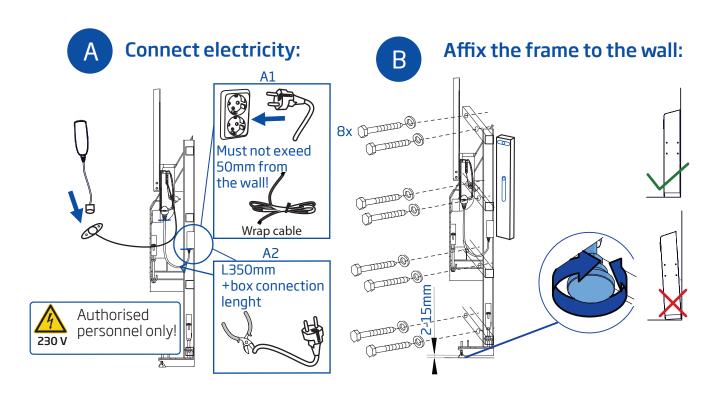


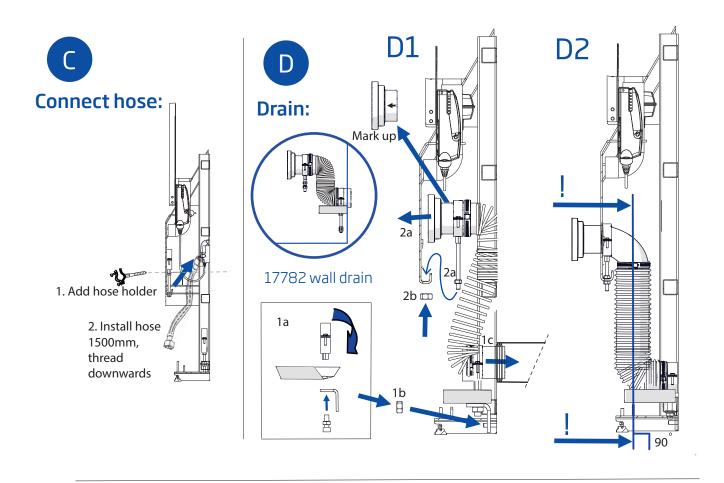


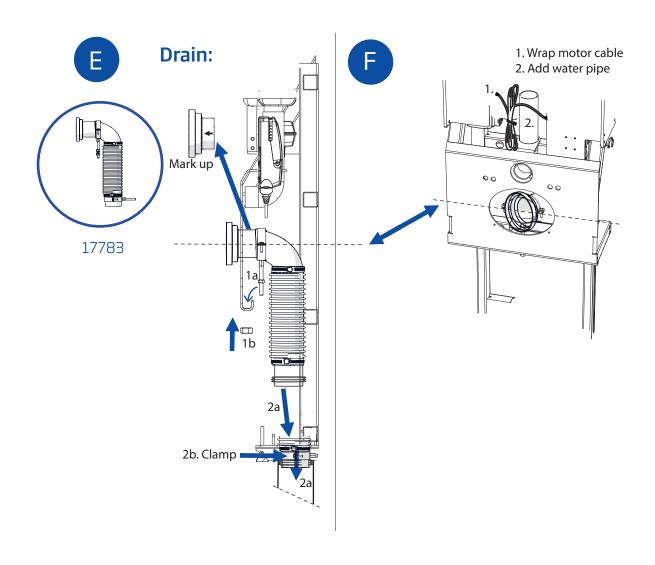


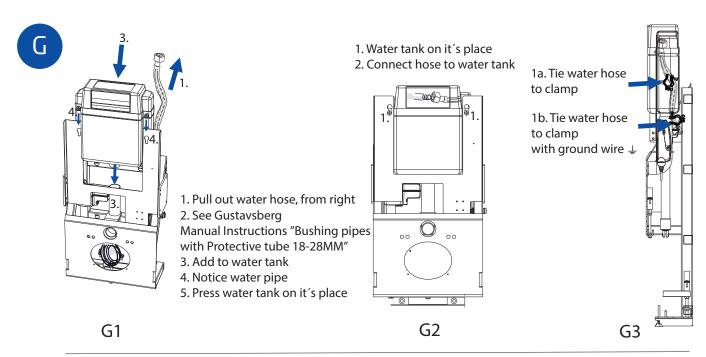


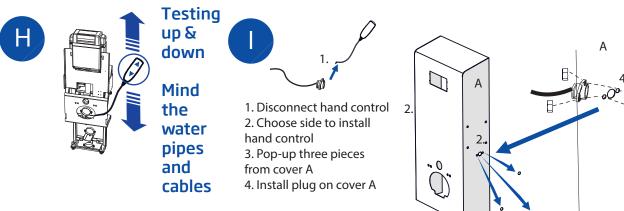

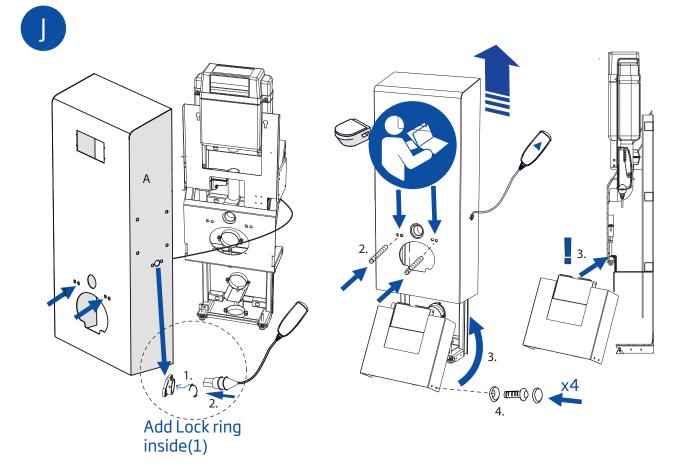


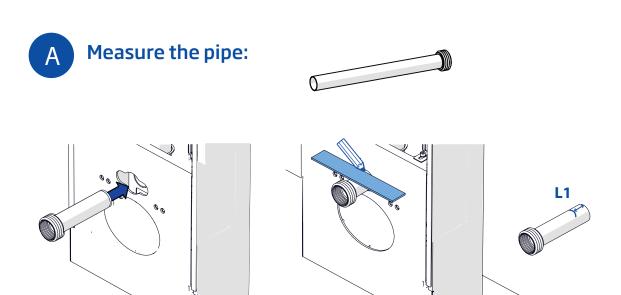


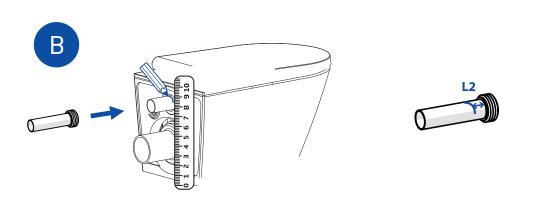

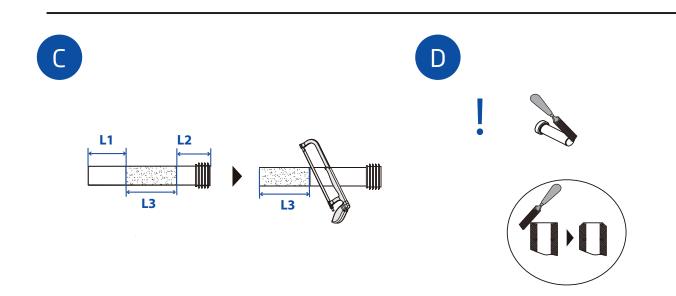



Mount the cover, seat and support rails:

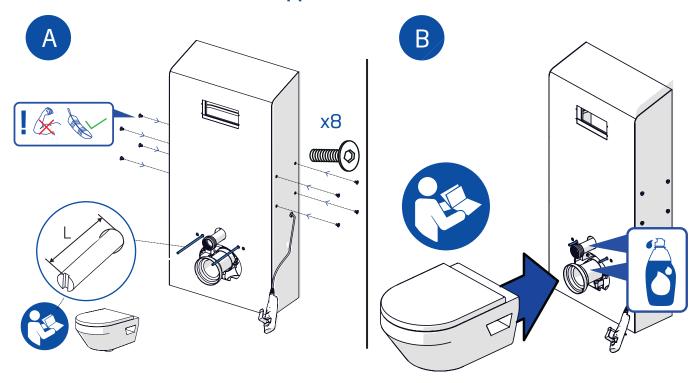

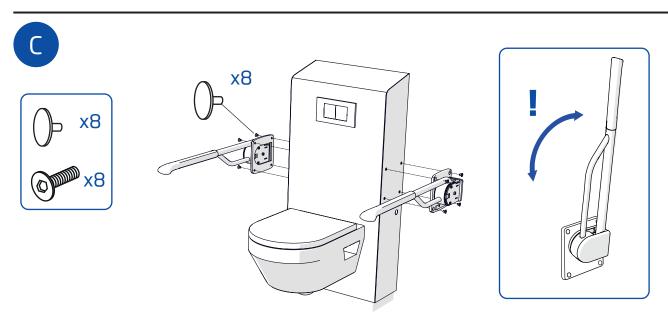

Mount the toilet roll holder, remote control and push button:

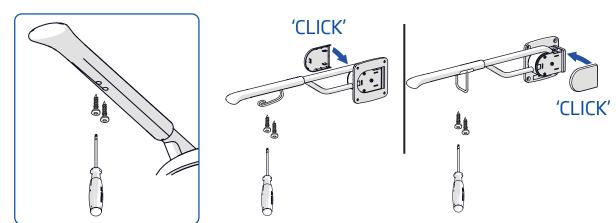


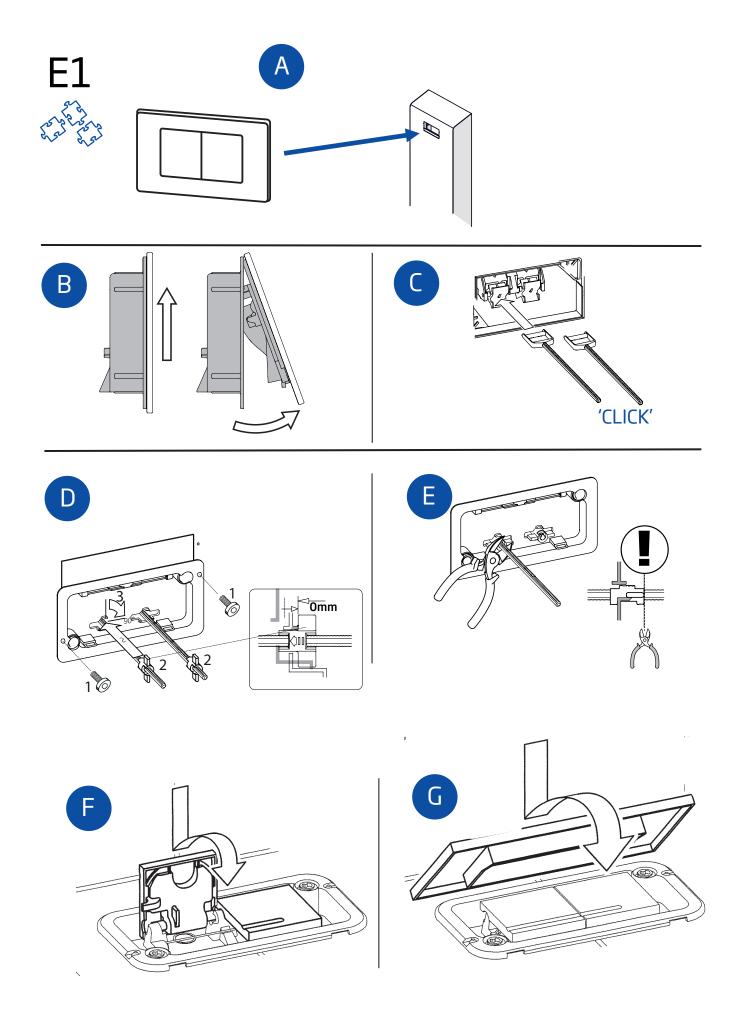





Always perform unclogging by moving the implement in the direction of the building's drain pipes.

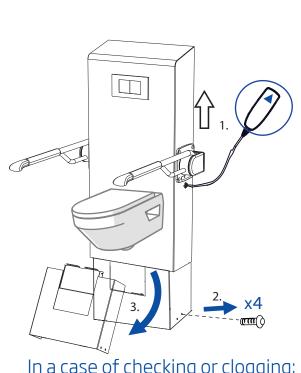

In a case of checking or clogging:

- 1. Lift Duo as up you can
- 2. Remove screws, 4 pcs
- 3. Remove under cover
- 4. Squeeze black drain pipe

